



Buoyancy Balance 6 kg X 0.1 g supplied with Support Frame Water Tank & Suspension Hook

Code: [42-1000/01](#)

Product Group: [Particle Density \(specific gravity\) and Water Absorption for Coarse Aggregates](#)

The buoyancy balance system developed by ELE consists of a rigid support frame, incorporating a water tank mounted on a platform. A mechanical lifting device is used to raise the water tank through the frame height immersing the specimen suspended below the balance. The balance supplied may also be used as a standard weighing device, thus providing a versatile and comprehensive weighing system in the laboratory. 6000 g x 0.1 g supplied with frame, water tank and suspension hook.

Standards

EN 1097-6, EN 12697-6, BS 812, ASTM C127, AASHTO T85

Further Information

****NOT AVAILABLE FOR SALE IN THE USA****


For 220 - 240 V AC, 50 - 60 Hz, 1 ph.

Accessories



200 mm Diameter Lid

Code: [79-0010](#)

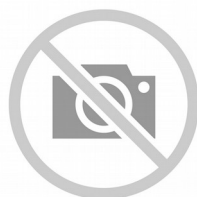

200 mm Diameter Receiver

Code: [79-0020](#)



200mm Dia BS Sieve 31.5mm Perforated Plate.

Code: 79-1630



200mm Dia BS Sieve 4.00mm Stainless Steel Mesh.

Code: 79-0310



200mm Dia BS Sieve 63 Mic Stainless Steel Mesh.

Code: 79-0070



200mm Dia BS Sieve 63.0mm Perforated Plate.

Code: 79-1670



Digital Pocket Thermometer

Code: [82-5420](#)



Sample Container - Tinned Steel 10 litres

Code: [81-3060](#)



Sample Tray - 406 x 406 x 50 mm

Code: [81-4030](#)



Wire Basket Brass with Handle Nominal 6000Cm³ Capacity with 1.7mm Wire Mesh.

Code: [42-1005](#)

Alternatives



ELE Buoyancy Balance

Code: [42-1000/02](#)